


# Ambu® Cardiac Care Trainer System W

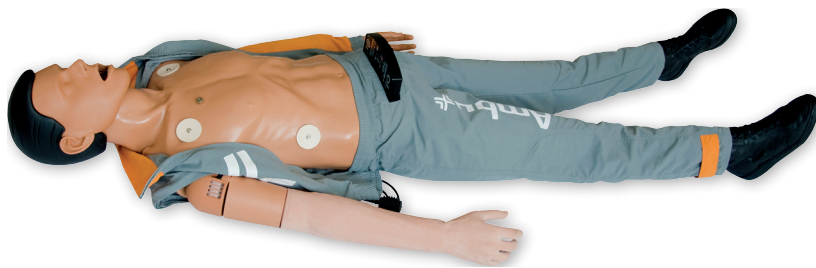

## Key Benefits

- Wireless connections
- Modular expandable
- The rigidity of the chest is adjustable
- Flexible electrodes
- Can be intubated
- PC-compatible
- ECG box is part of the software

## Ambu® Cardiac Care Trainer System W

The Ambu Cardiac Care Trainer System W (Wireless) is an advanced instruction- and training manikin that has been designed for training Advanced Life Support in ACLS (Advance Cardiac Life Support) algorithms and teamwork, either indoors or on location in field situations.

### Connection free

The Ambu Cardiac Care W does not need cables or wires, which means it can be placed away from computers and power supplies. This gives a great amount of freedom during training, and because there are no wires or cables lying around it gives an even more realistic scenario.

### Modular Expandable

Ambu Cardiac Care W is modular expandable, meaning that you can invest gradually and update the system as you go along.

### Mechanical Instrument

The Cardiac Care W is equipped with a mechanical monitoring instrument as well as digital sensors showing effectiveness of the resuscitation. It also shows results achieved for insufflations volume and depth of external chest compression, as well as it indicates any stomach inflation and incorrect hand positioning.

The Cardiac Care W can be connected to a computer through USB or WLAN using the Ambu CPR-Software for a more detailed analysis of the performed resuscitation.

The CPR Software is part of the wireless system in the Ambu Cardiac Care W. The CPR Software analyses and displays the performance done on the manikin in real time. In addition the CPR Software can display and control the ECG. The virtual ECG provides 26 rhythms and 2 artefacts and can automatically detect defibrillation.

# Specifications


Product	
Weight: Full body	13.3 kg
Torso with carrying case	21.1 kg
Full body with carrying case	25.1 kg
Full body incl. IV Trainer	25.9 kg
Total length: Full body	165 cm
ECG Box: length x width x height	18 x 14 x 4 cm
Weight: with batteries	545 g
Weight: with batteries	380 g
9V DC adapter: Length x width x height	8.6 x 5 x 7.2 cm
Weight	165 g




Ambu Cardiac Care in suitcase

# Ordering specifications

Item no.	Description
A284 317 000	Ambu Cardiac Care System W, Full body and Hard case with CPR-Software
A284 307 000	Ambu Cardiac Care System W, Full body, IV Arm and Hard case with CPR-Software
<b>Accessories</b>	
284 000 730	Ambu ECG Box
234 000 790	Track suit jacket
234 000 792	Track suit trousers
234 000 786	Legs with trousers and bag

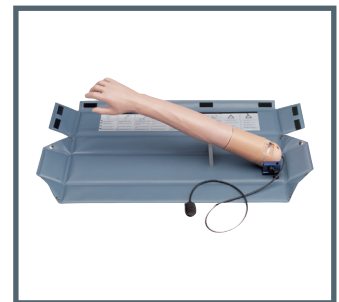

Ambu IV Trainer

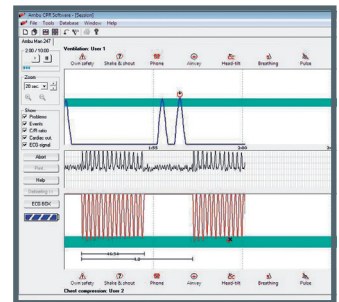

CPR Software screen